



Dobson Wedge Removal

After Removing the Lid...



Back of Wedge



Front of Wedge

Orient yourself with the front and the back of the Wedge

Remove the Three Set Screws



Remove the two set screws
on the front of the Wedge



Remove the one set screw at the
back of the Wedge

Remove Wedge Carefully




1. Grab the Wedge by the base



2. Gently lift up the back side of the Wedge with your other hand



3. Pull the Wedge through the front of the instrument

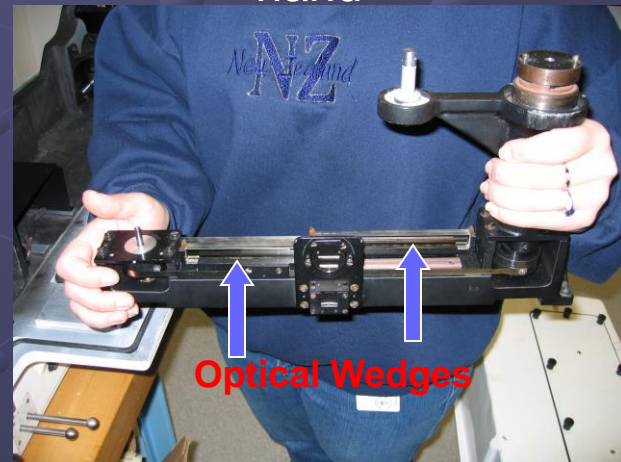


4. Hold the Wedge with both hands taking care not to touch the two optical wedges mounted within

Packing the Wedge

Avoid touching the optical bits, as oil from your skin can effect the calibration

Other optical Wedge

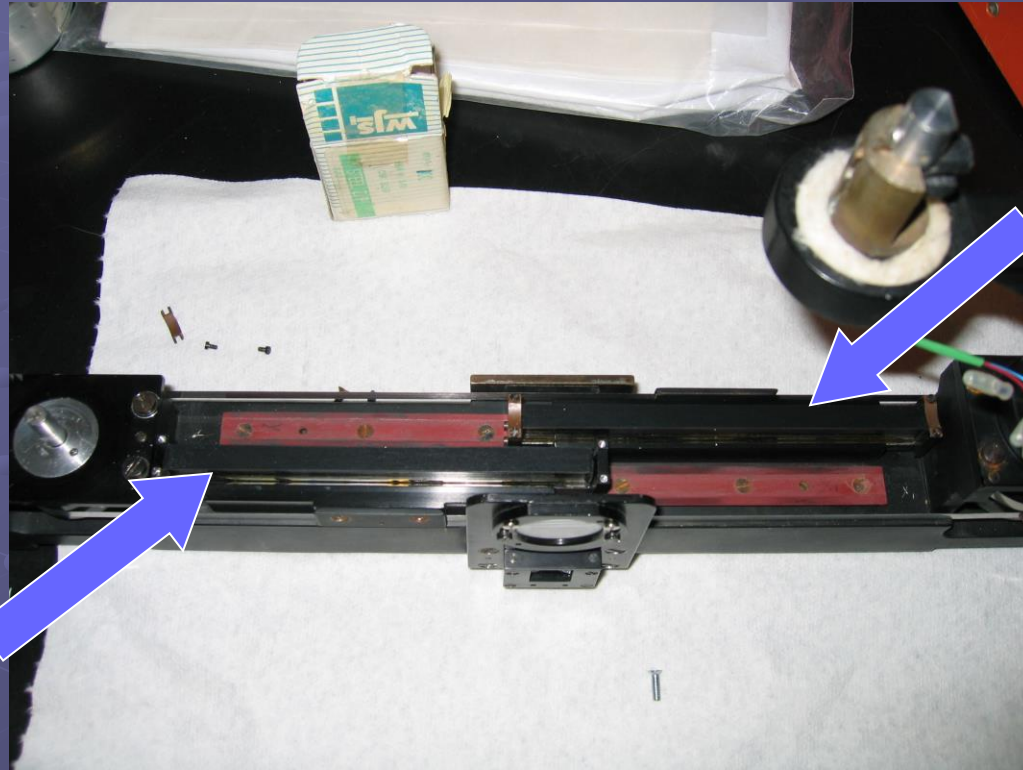

One optical Wedge

Take care that you do not touch the optical Wedges. The optical Wedges themselves need to be wrapped very carefully so that they are not damaged in shipping. They should be covered with several layers to protect them.

Help

● Help is here anytime you need it. Please call with any questions you may have!
Good Luck!!!

● Glen McConville 303-497-3989
720-282-9576

● Peter Effertz 720-314-8380